


Spore News

Process Challenge Devices: More Than a Biological Indicator in Fancy Packaging

Historically, medical device manufacturers and contract sterilizers had to place biological indicators (BIs) inside their medical devices or deep within a pallet of material in the most difficult to sterilize area so that they could adequately confirm the sterility of their products. Once the sterilization process was complete they had to retrieve the BIs and process them. There were numerous downsides to this process that made it an undesirable practice:

- Planting and retrieving the BIs was a time consuming process.
- The retrieval process compromised the sterility of the device or package it was planted in and increased the chance of contaminating other packages in the load.
- The technicians retrieving the BIs were potentially exposed to high levels of residual EO gas as they unpack the load to retrieve the BIs.

The best way to eliminate the need for planting BIs inside the load was to develop an external process challenge device (PCD) that could be placed on the outside of the load to mimic the performance of the internal BIs. Over the years, biological indicator manufacturers and end users alike developed PCDs to monitor specific loads. The down side of these devices was that they were a set design that provided a very specific resistance range. If the design didn't work for a particular load/process it was back to the drawing board. Mesa's line of PCD products with its broad range of available resistance levels and BI types eliminates this problem.

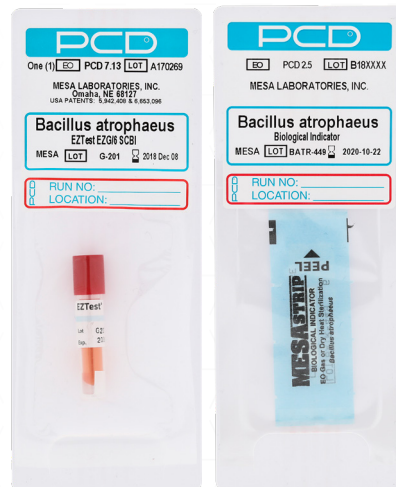

Mesa’s lineup of PCD products consists of standard biological indicators packaged in pouches composed of proprietary plastic film combinations. The pouch is attached to a plastic card that adheres to the exterior of the medical device packaging or pallet of material. Once a manufacturer or contract sterilizer qualifies the performance of a PCD configuration against their internal biological indicator performance and shows equivalence, they eliminate the need for the internal BI all together. All they have to do is affix the PCD to the exterior of the packaging/pallet. No more dismantling loads and compromising packaging while trying to locate and retrieve the BIs. Affixing the PCD to the exterior of the pallets allows for quick retrieval of the PCD, which significantly reduces the potential for EO exposure by the technicians, thus creating a much safer work environment.

Mesa currently offers six (6) PCD pouch types for monitoring ethylene oxide sterilization process. Each pouch type has a different permeability rate for ethylene oxide, thus allowing us to offer a broad range of resistances for a given type of biological indicator. The table below lists them in order of resistance from least to most resistant:

Table 1: Pouch Resistance Matrix

PCD Pouch	Resistance
Type 3.XX	Least
Type 5.XX	
Type 7.XX	
Type 6.XX	
Type 4.XX	
Type 2.XX	Most

Over the years you’ve heard us refer to biological indicators as a system. It’s not just the bacteria that creates the resistance of the product. It’s the bacteria, the inoculated substrate and the primary packaging that all come together to produce the D-value that we label all of our BI products with. If we apply the same “biological indicator is a system” approach to a pallet of product being sterilized with a BI placed at the most difficult to sterilize location, then we need to determine the resistance of this new system. As you could imagine, placing the BI in the interior of a pallet of goods significantly increases the resistance of a standard biological indicator. There are many different types of BIs available for monitoring EO sterilization, but even though there is a wide range of products available they all fall within a specific resistance range. That range is typically from 2.5-5.8 minutes with the overwhelming majority of commercially available product being in the range of 2.5-3.5 minutes. Once we take the packaging configuration of all the items on the pallet it’s possible that the resistance increases 10-fold.

This is because the sterilant has to pass through all of the packaging material before it gets to the BI. So in this instance, if an end user wanted to convert to an external PCD they would need something that produced a resistance of around 30minutes. Using our PCD products we are able to take a BI with a 3.0 minute D-value and offer a resistance range from as low as 3 minutes to as high as 58 minutes. By having such a broad resistance range we are able to provide the end user with a resistance that takes this new biological indicator into account. This is accomplished by selecting plastic films for the pouches that control the diffusion rate of EO through the PCD film.

The table below displays the nominal resistance for three types of common biological indicators packaged in the six (6) PCD pouch types:

Table 2: Resistance of Commonly Used BIs

PCD Type	BI Type	D-value
3.2	Naked Spore Strip	3min
5.2	Naked Spore Strip	4min
7.2	Naked Spore Strip	9min
6.2	Naked Spore Strip	12min
4.2	Naked Spore Strip	21min
2.2	Naked Spore Strip	28min

PCD Type	BI Type	D-value
3.5	Spore Strip in Glassine	4min
5.5	Spore Strip in Glassine	5min
7.5	Spore Strip in Glassine	14min
6.5	Spore Strip in Glassine	17min
4.5	Spore Strip in Glassine	32min
2.5	Spore Strip in Glassine	37min

PCD Type	BI Type	D-value
3.13	EZTest SCBI	3min
5.13	EZTest SCBI	13min
7.13	EZTest SCBI	26min
6.13	EZTest SCBI	35min
4.13	EZTet SCBI	50min
2.13	EZTest SCBI	58min

The tables listed on page 2 detail a number of different types of PCD: PCD type 2.13, 6.5, 3.2. What does that mean? All of the PCD products that we offer follow a unified catalog numbering format. The number starts with PCD and then is followed by a 2-3 digit number. The first number represents the type of pouch. If you order PCD 2.13 it will be composed of a type 2 pouch (represented by the first number). The 1-2 digit number following the decimal point represents the type of BI in the pouch. Below is an example key that shows the most common BI types used in PCDs:

Table 3: PCD Catalog # Format

PCD X.##	BI Description
1	Mesa BI Disc, 1-10006mm
2	Steris Log6 Spore Strip NA005
5	Mesa Log6 Spore Strip in glassine SGMG/6
6 or 62	3M Attest Log6 1264-S SCBI
7	Steris Log6 Spore Strip in glassine NA005
9 or 92	3M Attest Log6 1294-S Rapid Read SCBI
13	Mesa EZTest Log6 EZG/6 SCBI

Written by Robert Bradley, Sr. Director - Bozeman Operations, Mesa Labs

For additional questions on PCD selection or use,
contact bi-support@mesalabs.com

406.585.9535 or 877.746.8883 • biologicalindicators.mesalabs.com

Sterilization and Disinfection Control Manufacturing: 625 Zoot Way, Bozeman, MT 59718
Corporate Headquarters: 12100 West 6th Avenue Lakewood, CO 80228
ph: 800.992.6372 • 303.987.8000 • fx: 303.987.8989